



## PCO T-09 (Traje de Lycra sin mangas)


TALLE	CIRCUNSFERENCIA DE TRONCO (15)	AXILA - TROCANTE (IV)
1	48	30
2	54	34
3	54	38
4	57	42
5	60	46

Custom measurement form:  
**PCO ORTHOSES CARD**

PATIENT ID \_\_\_\_\_

PRODUCT CODE \_\_\_\_\_

NOTES \_\_\_\_\_

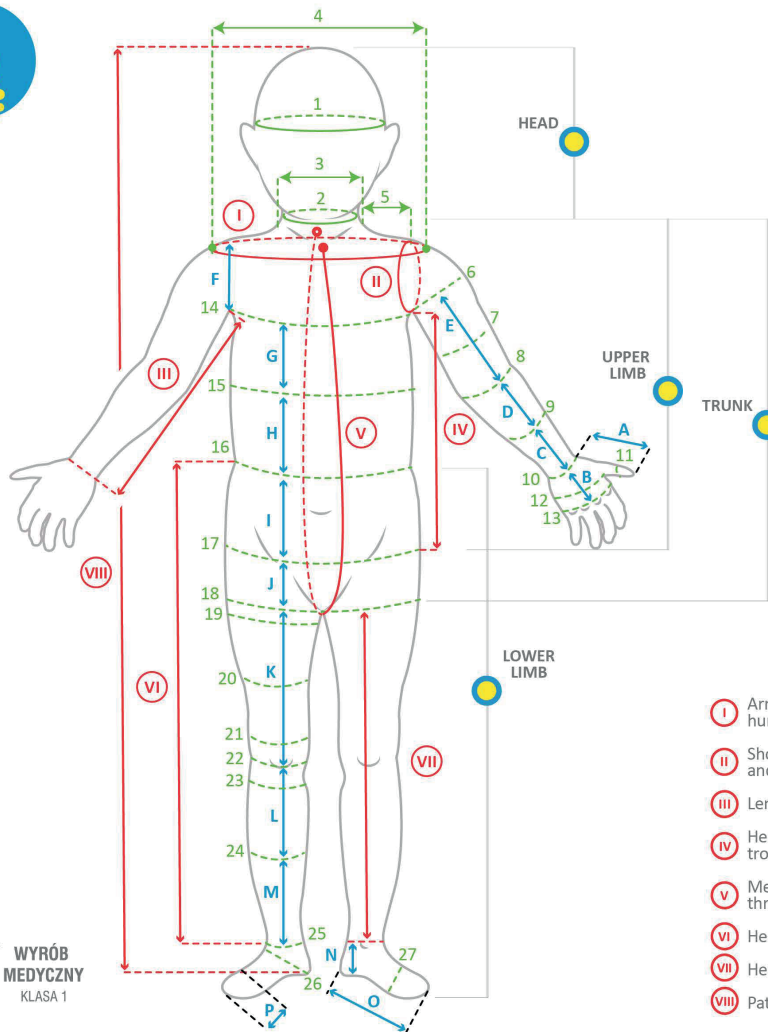


**⚠ WARNING!**

- Chest and upper limb measurements should be taken with patient lying down.
- Hip and lower limb measurements are taken with patient standing.
- All measurements should be taken in centimeters!
- If the finger put between the child's body and the PCO product meets resistance, the size of the product has been selected correctly.

**IMPORTANT!**

Only a correctly filled in measurement card is a guarantee of proper compression of the product. We are not responsible for wrongly taken measurements!

**HOW TO MEASURE?**

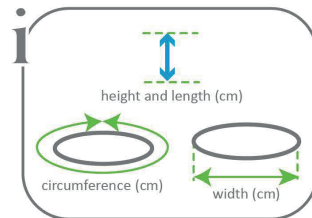

- I** Arms circumference at the level of the head of the humerus
- II** Shoulder circumference measured through the armpit and acromion
- III** Length measured from armpit to wrist
- IV** Height measured from the armpit to the greater trochanter of the femur
- V** Measurement from the manubrium of the sternum through the perineum to the cervical vertebra C7
- VI** Height measured from waist to ankle
- VII** Height measured from parineum to ankle
- VIII** Patient is height